
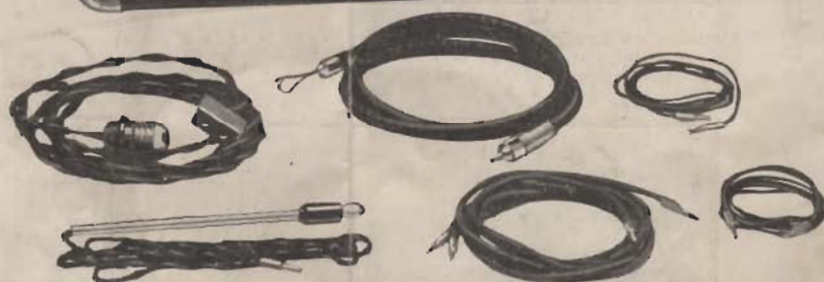


THE BETZCO ELECTRIC

A Superior Electro-Therapeutic Cabinet

Sold on Easiest Terms

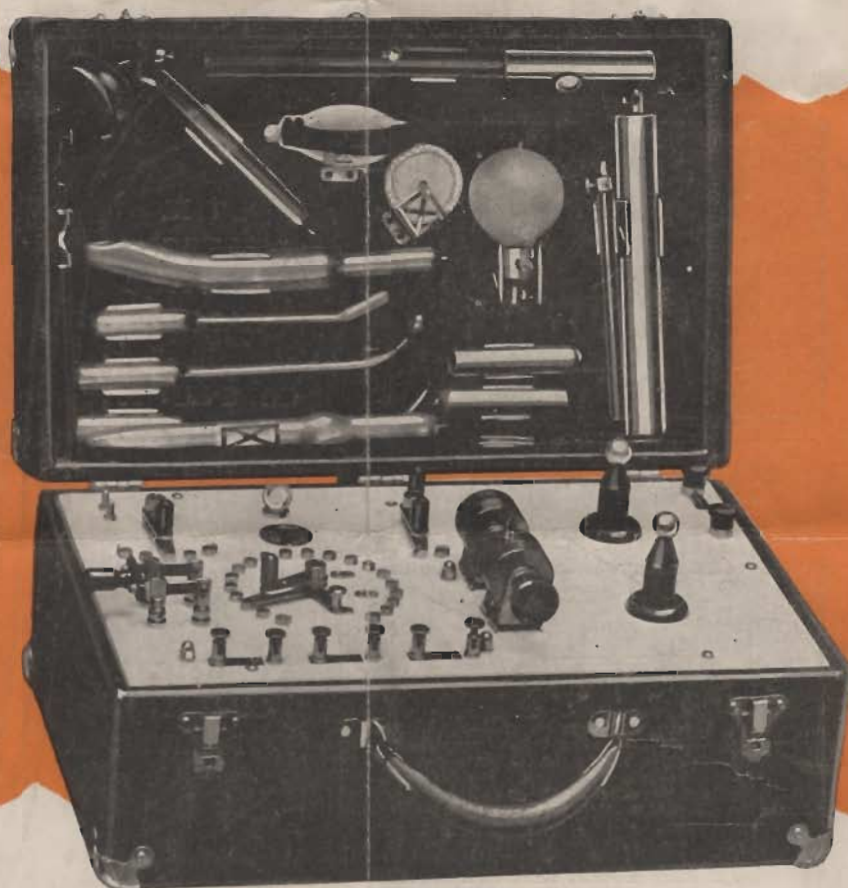

FRANKS. BETZ CO.

NEW YORK
6-8 West 48th St.

HAMMOND INDIANA

CHICAGO
30 E. Randolph St.

A COMPLETE AND COMPACT



Oral Cavity Examination

The Diagnostic Illuminator is equally valuable for all Orificial work. Connections 6 and 7.



Minor Nasal Surgery

The Cautery Handle and Knives furnished will be found a valuable aid for Sealing purposes, stopping Nasal Hemorrhages and opening Boils. Connections 8 and 9.



Sinusoidal Treatment

Patient's foot in tray of salt water; Spongio Disc just above knee, for testing Muscle Reflex. Connections 4 and 5.




Application Sinusoidal

One electrode on Lumbar Spine region—Spongio Disc between shoulders. Connections 4 and 5.

THE
BETZCO
ELECTRIC
CABINET
IS
COMPLETELY
EQUIPPED

OFFERS
ALL
PRACTICAL
CURRENT
MODALITIES

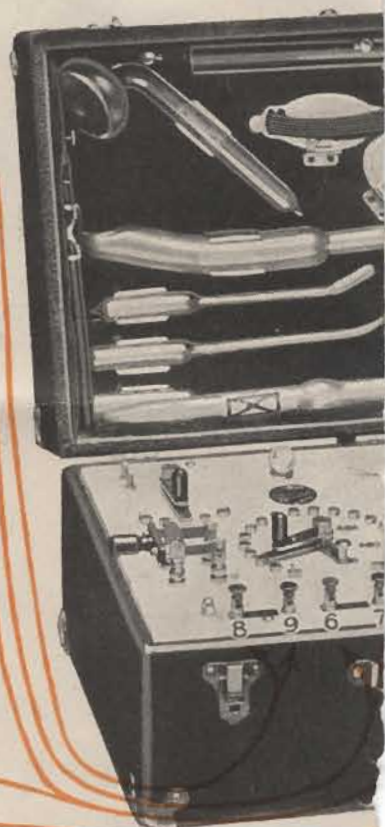

Vaginal High Frequency Treatment

Uni-polar treatment using Insulated Glass Electrode, with moderate current. Electrode in vagina. Connection 3.



Rectal High Frequency Treatment

Insulated electrode in rectum, using uni-polar method of application. Very valuable for relieving Prostatic trouble. Connection 3.

HIGH FREQUENCY CABINET




Abdominal Diathermy Application

Using the 8-inch Flexible Pad on the Abdomen and a large Block Tin electrode under the Spine. Connections 1 and 2.



A "Modified Diathermy" Application

The operator uses his fingers as a Rheostat. Treatment relieves congestion, and produces an active Hyperemia wherever used. Connections 1 and 2.

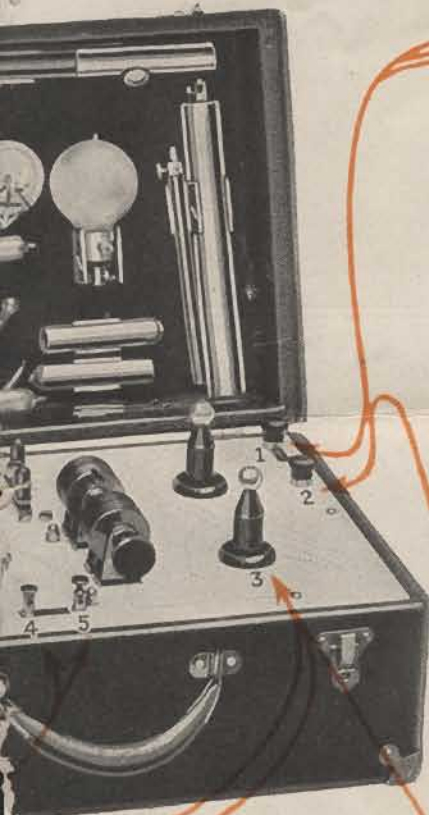

A
SUCCESSFUL
OUTFIT
FOR
GENERAL
OFFICE
PRACTICE

REMARKABLY
EASY AND
SIMPLE
FOR THE
OPERATOR



Diathermy Application to Chest

Large indifferent electrode under spine and flexible sheet of Block Tin (6-inch diameter) on chest. Connections 1 and 2.



Heating the Inner Ear

An application of Diathermy Current—electrodes in patient's hands, with index fingers in ears, causing heating through the head. Connections 1 and 2.




Uni-Polar Fulguration

The Fulguration Electrode with the curved tip connected to the High Frequency outlet. Mild current. Very effective for the removal of Warts, Moles, Lupus Vulgaris, etc. Connection 3.



Auto-Condensation

Patient on Chair Pad, holding the metal handle in both hands. Especially valuable for the reducing of High Blood Pressure. Connections 1 and 2.

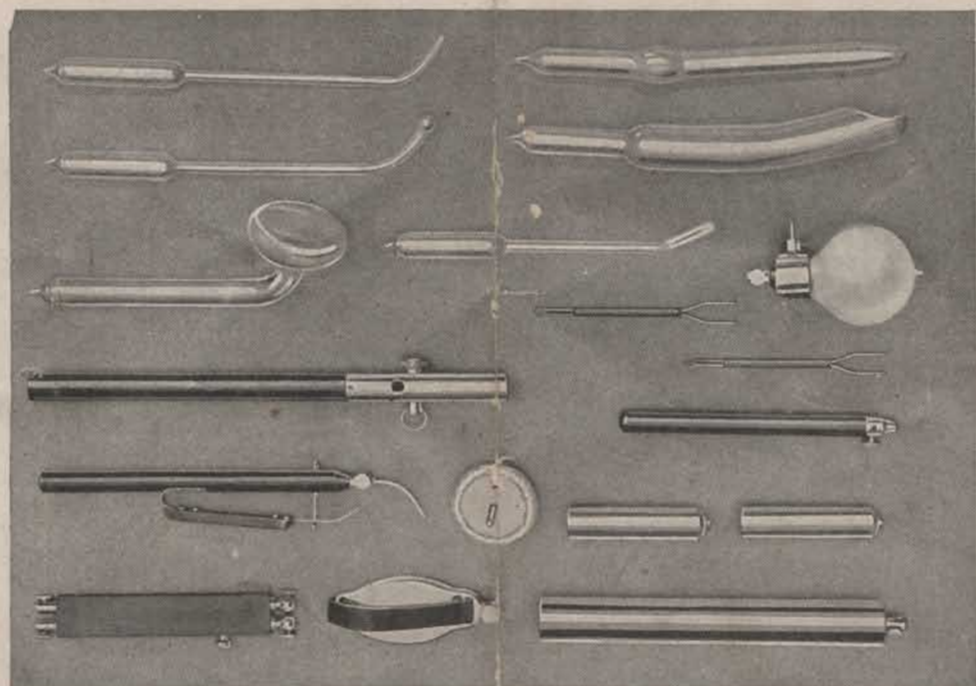

Completely Equipped with Modern Accessories

The Betzco Electric is a powerful, simple, practical and durable cabinet designed to fill the needs of general office practice. The use of this outfit gives you an opportunity to handle cases that rightfully belong in your office, but are now being handled outside of the profession. Cabinet offers the following modalities:

High Frequency—Its effects are noted in increased blood supply, increased oxidation, secretions and local nutrition.

D'Arsonval—Indicated in treatment of high blood pressure and for elimination through cellular massage.

Auto Condensation—Delivered through chair pad included with outfit. A soothing treatment for the entire system.

Diathermy—A high frequency current which produces heat effects deep in the tissues. Heat may be directed to the exact area desired and intensity of application modified to suit the case.

Fulguration—For desiccating or cauterizing in the removal and treatment of warts, moles, growths, including destruction of papilloma of the bladder.

Modified Diathermy—Used for relieving congestion and for producing an active hyperemia.

Ozone—Indicated in treatment of infectious diseases, anemia, diseases of the respiratory organs and in conditions where there is imperfect oxidation and impaired nutrition. Special electrodes for inhalation can be furnished separately.

Sinusoidal—Indicated in the treatment of paralysis and habitual constipation. Sometimes used effectively in treatment of deafness.

Cautery—Sufficient cautery current generated to heat any of the ordinary cautery knives. Has 30 points of control.

Diagnostic—Delivers an accurately regulated current for lighting diagnostic lamps from 1/2 to 12 volts. Operates the smallest cystoscopic lamps safely.

The following accessories included with the Betzco Electric:

One Latest Model Generator in leatherette covered case, velvet lined with marble top and nicked clips for equipment.

One Diathermy Hand Electrode.

One Connecting Cord and Plug.

One Book of Complete Directions.

One Adjustable Electrode Handle.

One Surface High Frequency Electrode.

One Rectal High Frequency Electrode.

One Vaginal High Frequency Electrode.

One Throat High Frequency Electrode.

One Nasal High Frequency Electrode.

One Urethral High Frequency Electrode.

One Heavy Insulated Cord.

One Auto Condensation Pad.

One Regulating Lamp and Clamp.

One Milliampere Color Scale.

One D'Arsonval Handle.

Two Sinusoidal Cords.

Two Sinusoidal Handles.

One Diagnostic Illuminator with cords, holder and Tungsten lamp.

One Cautery Handle, adjustable.

One Cautery Knife, pointed.

One Cautery Knife, looped.

Two Cautery Cords.

One Fulguration Electrode.

One straight Tip for Fulguration Electrode.

One Curved Tip for Fulguration Electrode.

One Flat Tip for Fulguration Electrode.

One Felt Sinusoidal Disc and Handle.

ORDER COUPON

FRANK S. BETZ CO.,

HAMMOND, INDIANA.,

DEAR SIR:

Enclosed is remittance for \$20.00 for which you may send me the 9A213 Betzco Electric High Frequency Cabinet. I can return it for full credit if not satisfactory in every respect.

NAME

ADDRESS

CITY STATE

**YOURS ON
EASIEST TERMS**

Send \$20.00 Only

The 9A213 Betzco Electric Cabinet regularly sells for \$110.00. We offer you special terms, however, of \$20.00 only with your order and the balance of \$90.00 in equal monthly payments of \$10.00 each without interest. Extra fees collected for treatments and applications should easily pay for the outfit.

WITH INSTRUCTION BOOK

Outfit includes a complete and detailed book of instructions from which you can quickly learn to operate this outfit effectively. This booklet is written in simple terms and is very easy to interpret.